INTEGO.TRANSCENDENTAL

TREATMENT CENTER





RISE ABOVE THE ORDINARY WITH INTEGO

Based on Sirona's heritage of combining innovative technology and award-winning design, INTEGO is an entirely new generation of treatment centers that is as admirable as it is affordable. Reliable and efficient stylish and contemporary, INTEGO proves that excellence and value are not mutually exclusive.



TRANSCENDENTAL

A treatment center unlike any other, INTEGO takes the entire treatment experience to new heights of success through the thoughtful synthesis of innovation, style, and technology. INTEGO represents dentistry at its most extraordinary, distinguishing your practice as a pioneer among your peers.

INTEGRATION



FLEXIBLE CONFIGURATIONS FOR THE WAY YOU WORK.

Your workday is complex and places many different demands on the treatment center you choose. INTEGO provides the perfect solution to integrate with your workflow: A choice between two models with flexible configurations designed to meet your needs, each offering optimal performance.



INTEGO TS
Traditional Delivery



INTEGO CS
Continental Delivery

BOTH MODELS OFFER THE FOLLOWING:

- EasyTouch user interface, allowing complete control of the instrument in use, with a screen showing you what you need when you need it.
- The USB port offers you the possibility to connect your Schick sensor and take x-ray images right at the treatment center with the optional HELIODENT^{PLUS} intraoral source.
- Integration with SIDEXIS 4 with direct control through the EasuTouch touch screen
- The Endodontic function allows you to use your file system with the INTEGO's interface EasyTouch, and the ENDO 6:1L attachment.
- 22" HD monitor that allows you to explain intraoral and x-ray images, software and planning views, media player files, and PowerPoint presentations to your patients directly at the treatment center.
- A network interface in the patient chair and a USB port in the dentist or assistant element allow
 you to integrate new functions and network your practice—troubleshooting by your service
 technician is also much easier



I DESIGN

CREATED TO BE WITHIN YOUR REACH

Brilliance is the fusion of style and simplicity, the combination of contemporary design components at work, and the connection between you and your patients. From the thoughtful ergonomic features, to the reliable armrests, to the almost visceral elevation system—INTEGO provides optimal comfort to you and your patients throughout the treatment process.

HIGH QUALITY

Durable materials, quality workmanship, and sturdy, future-proof design.

CLEAN LINES

Flowing, continuous surfaces that are as pleasing to the eye as they are easy to clean.

THIN BACKREST

Allows easy and unencumbered access to the oral cavity.

UNENCUMBERED MOVEMENT

Smooth chair back and seat bottom; even the right armrest can be rotated forward for better access.





EXPERIENCE

PATIENTS COME IN ALL SHAPES AND SIZES, BUT COMFORT IS UNIVERSAL.

The INTEGO patient chair features an ergonomic design and, together with the headrest, ensures that your patients are especially comfortable. The chair lowers far enough for even shorter patients to get in and out easily. The double articulating headrest is quickly adjustable to various treatment situations and patients with different statures, and the variable foot rest allows taller patients to rest their feet comfortably. High quality upholstery enhances comfort even during extended treatments, and the tilting backrest matches the movement of the patient's spine to prevent stretching and compression effects. One thing is certain—patients who realize a comfortable and relaxed experience in the chair feel more secure and remain loyal to your practice.



INTEGO EQUIPMENT OPTIONS

OPTIONS FOR THE PATIENT CHAIR

2 treatment positions, S and O position	
Number of user-specific settings	2
Chair lift system 407 lbs (185 kg)	
Double articulating headrest	
Premium upholstery	e de la
Lounge upholstery	
4-way foot switch	
Booster cushion for small patients	
Armrests	
Pressure reduction kit	
Cover installation kit	
Separate junction box	
Base plate for uneven floors	

WATER UNIT

Cuspidor with manual swivel function	
Fresh water bottle with switch to main water supply	
Water flush hose cleaning system	

DENTIST ELEMENT

Height-adjustable support arm with pneumatic brake	
5 instrument positions	
3 Drill drives: Turbines or up to 2 electric	
Pneumatic instrument foot control	
Advanced electronic foot control	
EasyTouch user interface	
SPRAYVIT E with light (water on the left)	
Standard 3-way syringe	
Turbine hose with light	
BL ISO E micro-motor with hose	
SIROSONIC L scaler	
Intra-oral light source kit for turbine and SPRAYVIT E	
Endo function	
Purge and manual sanitization function	
CS tray holder for 1 standard size tray	
Additional TS tray holder	
Removable silicone sleeves for handles	

ASSISTANT ELEMENT

Comfort version, on a support arm with wide swivel range with 4 instrument positions	
Saliva ejector and HVE	
SPRAYVIT E with light (water on the left)	
Standard 3-way syringe	
2nd HVE kit	
Touch keys	
Silicone pad for Comfort assistant element	

OPERATING LIGHT OPTIONS

LEDlight	L
LEDview	
HELIODENTPLUS mounted on light post	

PATIENT COMMUNICATION

USB interface in the dentist element	
Monitor mount on the light support column	
Network interface	

As the heart of your practice, INTEGO harmoniously combines appearance and functionality. In addition to its harmonious design, coordinated components, and a variety of upholstery selections, INTEGO provides optimal access to the patient and smooth treatment workflows with a high level of patient comfort.

Standard



Optional



LEDVIEW OR LEDLIGHT OPERATING LIGHTS

Ideal illumination of the treatment area.




HELIODENTPLUS

Readily integrable with INTEGO, the HELIODENT^{PLUS} x-ray tube head fits into every practice and can be mounted on the wall, ceiling, or as a mobile unit.



DOUBLE ARTICULATING HEADREST

Easier to position and allows optimal visibility.

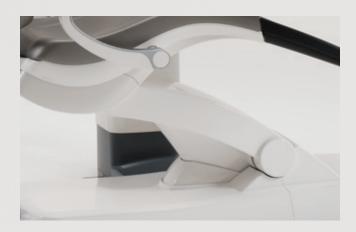

ELEVATION SYSTEM

Gentle chair lift drive accommodates patients up to 407 pounds.



CUSPIDOR

Removable and swiveling, easy to clean, and easily accessible.



4-WAY FOOT SWITCH

Regulates up/down and forward/backward movement of the chair, preset patient chair positions, and the advanced electronic foot control—the central control unit of INTEGO.

CALL YOUR PATTERSON REPRESENTATIVE FOR MORE INFORMATION.

Your Patterson Representative will help answer any questions and guide you through the process of getting started with INTEGO.

800.873.7683 SIRONATC.COM

Exclusive Distributor



SIRONA DENTAL, INC. 4835 Sirona Drive Charlotte, NC 28273

